

# VivaSight-DL

The new generation of double lumen tubes for Single Lung Ventilation

## Key Benefits

- Provides real-time visual monitoring during placement and throughout the entire SLV procedure
- Helps overcome today's placement and positioning challenges of SLV
- Makes lung isolation fast and effective
- Provides a clinical benefit and improves patient safety
- Integrated with Ambu® aView™ monitor



## VivaSight-DL

### Continuous monitoring during Single Lung Ventilation

Intubation and airway management during isolation procedures present a series of challenges for even the most experienced team. Visualisation during Single Lung Ventilation (SLV) increases patient safety and is increasingly recommended as good clinical practice.

### Single Lung Ventilation made simple and safe

VivaSight-DL is the new generation of Double Lumen Tubes with integrated high-resolution camera providing placement of the DLT fast and effective as well as continuous visual

monitoring throughout the procedure. Malpositioning and dislocation are easily detected and corrected. The real-time high resolution video image is transmitted to Ambu® aView™ monitor. Secretion or blood can be cleaned using the integrated flushing system.

The full Ambu portfolio for SLV also includes VivaSight-SL (ETT) with integrated high-resolution camera together with the bronchial blocker VivaSight-EB, and the single-use bronchoscope aScope™ 4 Broncho Slim.

# Specifications

VivaSight-DL with aView adapter cable	35 Fr Left	37 Fr Left	39 Fr Left	41 Fr Left
Size	35 Fr	37 Fr	39 Fr	41 Fr
ID (mm- effective)	4.0	4.2	4.8	5.0
OD (mm- bronchial)	10.5	11.5	11.5	12.0
Length (±10 mm)	330	330	330	330


VivaSight-DL Set Includes	Quantity
VivaSight-DL with aView adapter cable	1
Suction Catheter Straight 8 FR	2
Suction Catheter Curved 10 FR	2

General
Usage: Single-use only
Pack: Supplied individually packaged, ETO-sterilized
Shelf life: 36 months
Storage: store at temperatures between 0°C- 42°C in a dry and dark place
Materials
Tube: Medical-grade PVC; additional injection port for rinsing camera optics
Cuff: Medical-grade PVC, high volume, low pressure
Tube hardness: 87 (Shore A)

Camera
Image sensor: CMOS
Resolution: CIF 320 x 240 (76,800 pixels)
Depth of field: 12-60 mm
Field of view: ~100° diagonal
Light Source
LEDs: 2 white
LED power: 49.5mW - 66mW
Electrical
Power: Up to 36 mA @ DC 5V

## Ordering specifications

VivaSight-DL with aView adapter cable	Order number	Quantity/box
VivaSight-DL 35 Fr - Left	DLVT35LAS	5
VivaSight-DL 37 Fr - Left	DLVT37LAS	5
VivaSight-DL 39 Fr - Left	DLVT39LAS	5
VivaSight-DL 41 Fr - Left	DLVT41LAS	5
Monitor	Order number	Quantity/box
aView monitor	405002000	1

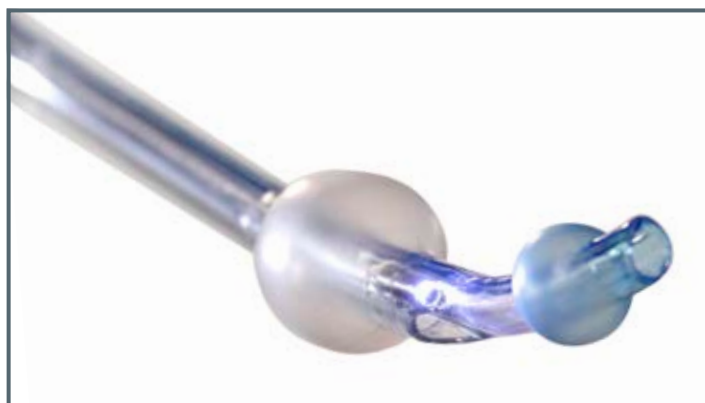

VivaSight-DL



VivaSight-DL connected to Ambu® aView™ monitor

Manufactured by **ETView Ltd.**  
 Catom 2 Street  
 Misgav Business Park  
 M.P. Misgav 2017900  
 Israel

VivaSight is distributed by **Ambu A/S**  
**Ambu A/S**  
 Baltorpbakken 13  
 DK-2750 Ballerup  
 Denmark