

VivaSight-EB



Key Benefits

- Compatible with VivaSight-SL for a complete airway management solution

VivaSight-EB

Single Lung Ventilation made simple and safe

VivaSight-EB is a single-use endobronchial blocker for Single Lung Ventilation (SLV). For visual guidance and positioning the VivaSight-EB can be used with VivaSight-SL single lumen tube with integrated high-resolution camera.

The VivaSight-EB consists of a sterile, single-use, “steerable” balloon-tipped catheter that is visually guided to a selected bronchial segment to affect balloon blockage of the distal airway.

The angled distal tip of VivaSight-EB can be adjusted to facilitate navigation to the bronchi. When used in conjunction with the VivaSight-SL, continuous monitoring throughout the procedure ensures that dislocation can be easily detected and corrected.


For Single Lung Ventilation the Ambu visualisation portfolio includes VivaSight-DL, a double lumen tube with integrated camera, VivaSight-SL, a single lumen tube with integrated camera, VivaSight-EB, as well as the single-use bronchoscope aScope™ 4 Broncho Slim.

Specifications

VivaSight-EB	Size	Length (±5 mm)
EBB09700AS	9Fr	695 mm
<p>General: VivaSight-EB features a blue silicone cuff.</p> <p>Usage: Single-use only. Sterile (ETO) unless package has been opened or damaged.</p> <p>Pack: Supplied individually packaged.</p> <p>Storage: Store in a dry, cool and dark place.</p> <p>Shelf life: 5 years.</p>		

Ordering specifications

VivaSight-EB	Order number	Quantity/box
VivaSight-EB	EBB09700AS	5
aScope 4 Broncho	Order number	Quantity/box
aScope 4 Broncho Slim 3.8/1.2 mm	476001000	5
VivaSight-SL with aView™ adapter cable	Order number	Quantity/box
VivaSight-SL 7.0 mm	TVT70100AS	5
VivaSight-SL 7.5 mm	TVT75105AS	5
VivaSight-SL 8.0 mm	TVT80110AS	5
VivaSight-SL with VivaSight-EB 9Fr Pack	Order number	Quantity/box
VivaSight-SL 7.0 mm + VivaSight-EB 9Fr	TVEB70100AS	5
VivaSight-SL 7.5 mm + VivaSight-EB 9Fr	TVEB75105AS	5
VivaSight-SL 8.0 mm + VivaSight-EB 9Fr	TVEB80110AS	5
Monitor	Order number	Quantity/box
aView	405002000	1

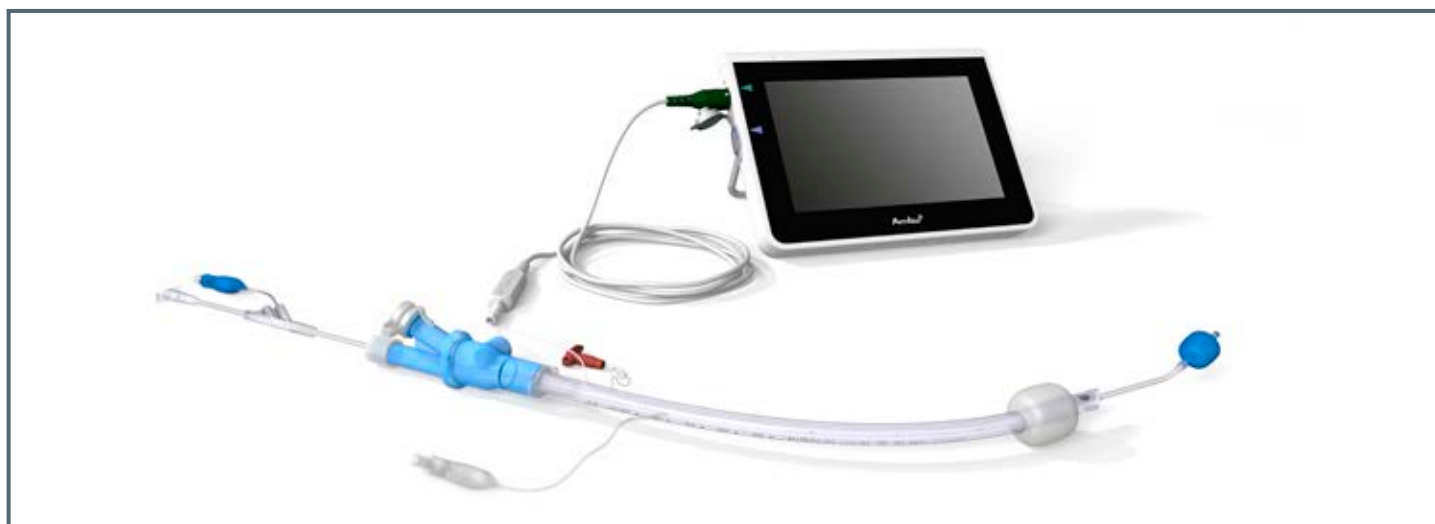

VivaSight-EB used with VivaSight-SL and aView monitor

Manufactured
by **ETView Ltd.**
Catom 2 Street
Misgav Business Park
M.P. Misgav 2017900
Israel

VivaSight is distributed by **Ambu A/S**
Ambu A/S
Baltorpbakken 13
DK-2750 Ballerup
Denmark

Rx ONLY **CE** 0483