

Ambu® aFlow

Pediatric Coaxial Circuits




Coaxial anaesthesia circuit are reliable convenient circuits designed for enhancing natural warming of inspired gases. The circuits have a streamlined tube-within-tube design, which reduces drag on the patient airway. The system provides increased circuit mobility and has a low weight compared with a traditional circle circuits. The coaxial system is ideal for use within both the operating theatre and ICU. It is well suited for low-flow anaesthesia.

Key Benefits

- Lower tube compliance
- Transparent tubing which makes viewing inside easier
- Smooth inner surface for easy drainage of condensed water
- Soft flexible connectors ensure a perfect seal with the anaesthesia machine inlet and outlet

Pediatric Coaxial Circuits

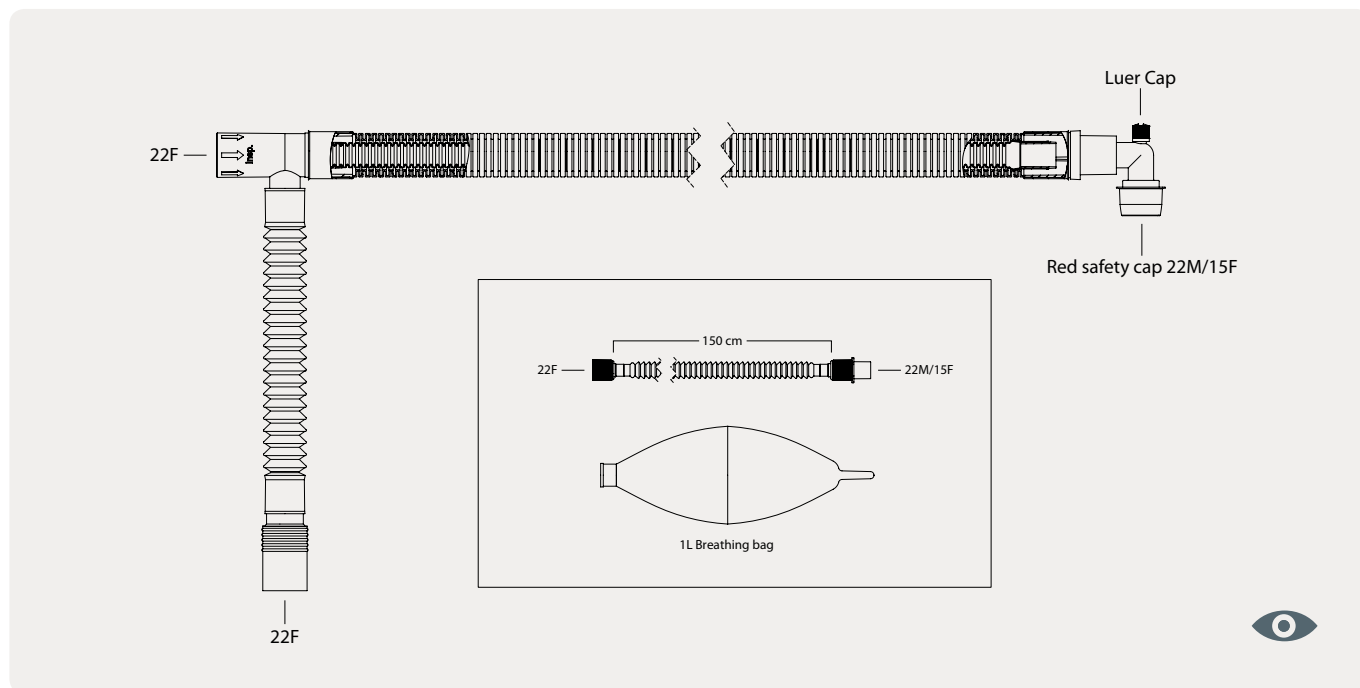

Ambu item no.	350208600	350208700
Size	Pediatric	Pediatric
Quantity/Box	25	15
Material	PE	PE
Tube length	1.6 m	1.6 m
Tube O.D/I.D.	22 mm / 15 mm	22 mm / 15 mm
Elbow Connector with luer cap	✓	✓
Red safety cap	✓	✓
Extention limb	-	1.5 m
Breathing bag, Volume	-	1 L
Compliance @60 cmH ₂ O	0.75	0.7
Resistance to Flow @30 L/min	0.045	0.045

Distributed by:
Ambu A/S
 Baltorpbakken 13
 DK - 2750 Ballerup
 Denmark
 Tel. +45 7225 2000
 Fax +45 7225 2053
 www.ambu.com

Manufactured by:
Altera Tibbi Malz.
San. ve Tic. A.S.
 İbni Melek OSB Mh. Tosbi Yol 4 Sok. No:29
 35900 Tire İzmir
 TURKEY
 Tel: (+90 232) 513 51 10 (Pbx)
 www.altera.com.tr