

Ambu[®] aFlow

Adult/Pediatric Septum Divided Circuits




Circuits that are suitable for use in both anaesthesia and ICU.
Available with and without water traps and in different lengths.
The circuit is lightweight and divided into two lines for expiration and inspiration, by a longitudinal flexible septum that runs through the tube.

Key Benefits

- Septum allows the transfer of heat from the expiratory to the inspiratory line
- Inspiratory air is naturally warmed and hence moistened

Adult/Pediatric Septum Divided Circuits

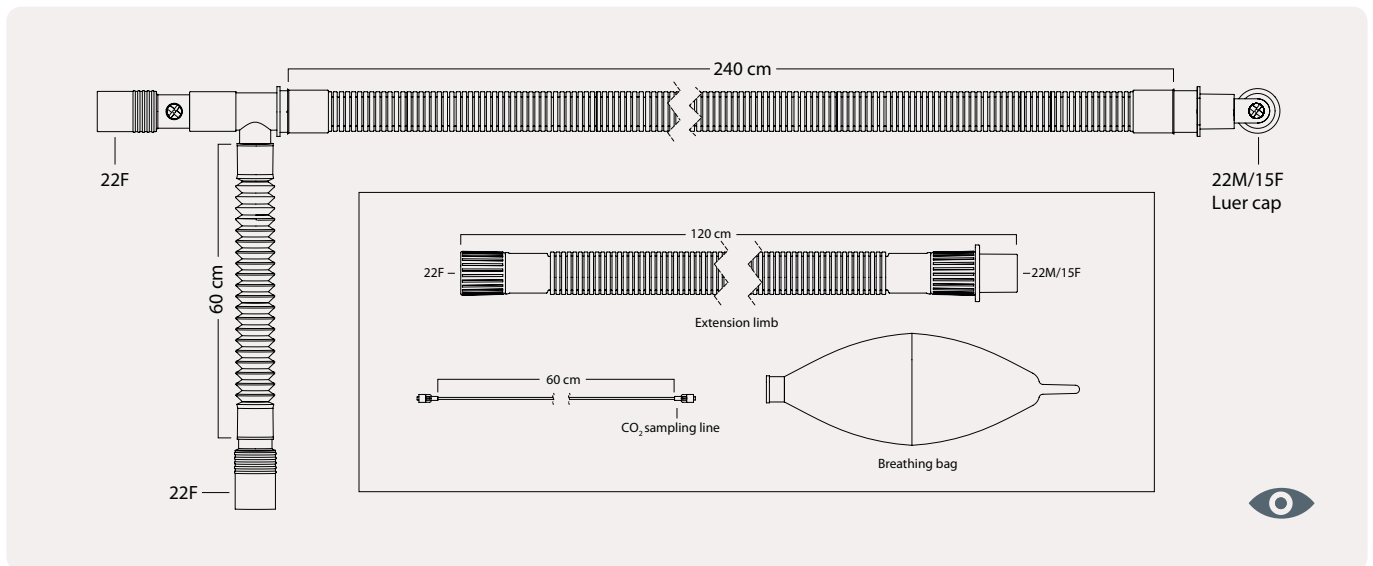

Ambu item no.	350201200	350201300	350201400	350201500	350201600	350201700	350201800
Size	Adult/ Pediatric	Adult/ Pediatric	Adult/ Pediatric	Adult/ Pediatric	Adult/ Pediatric	Adult/ Pediatric	Adult/ Pediatric
Quantity/Box	15	15	20	20	15	10	10
Material	PE	PE	PE	PE	PE	PE	PE
Tube length	1.8 m	2.4 m	1.8 m	3.0 m	1.8 m	1.8 m	1.8 m
Tube I.D.	22 mm	22 mm	22 mm	22 mm	22 mm	22 mm	22 mm
Elbow Connector with luer cap	✓	✓	✓	✓	✓	✓	✓
Double Swivel Elbow Connector Head	✓	✓	✓	✓	✓	✓	✓
Red safety cap	✓	✓	✓	✓	✓	✓	✓
Water trap	-	-	-	-	-	Integrated	Integrated
Extention limb	1.2 m	1.2 m	1.5 m	1.5 m	1.2 m	-	1.5 m
Breathing bag, Volume	2 L	2 L	-	-	2 L	-	2 L
CO ₂ Sampling line integrated	✓	✓	-	-	-	-	-
CO ₂ Sampling line	0.6 m	0.6 m	-	-	-	-	-
Temperature port, machine end	✓	✓	-	-	-	-	-
Compliance @60 cm H ₂ O	0.87	1.06	0.86	1.3	0.84	0.99	0.99
Resistance to Flow @30 L/min	0.071	0.0845	0.0605	0.079	0.0585	0.0545	0.0545

Distributed by:
Ambu A/S
 Baltorpbakken 13
 DK - 2750 Ballerup
 Denmark
 Tel. +45 7225 2000
 Fax +45 7225 2053
 www.ambu.com

Manufactured by:
Altera Tibbi Malz.
San. ve Tic. A.S.
 İbni Melek OSB Mh. Tosbi Yol 4 Sok. No:29
 35900 Tire İzmir
 TURKEY
 Tel: (+90 232) 513 51 10 (Pbx)
 www.altera.com.tr