

# Ambu® SureSight™ Connect

Maximise your focus on the patient

## DATASHEET

SureSight Connect is a video laryngoscope, which consolidates our trusted airway visualization platform. It can help you to maximise your focus on the patient as you move through the airway management algorithm.

### Designed to save space and reduce task complexity

The Ambu® SureSight™ Connect video laryngoscope solution puts your important, familiar tools within reach and ready to use. The platform is comprised of a single-use video laryngoscope, bronchoscope, rhinolaryngoscope and OLV solution all compatible with the same endoscopy system (displaying and processing unit).

### Use two devices on the same screen simultaneously

With the dual view feature in aView™ 2 Advance or aBox™ 2\*, you are able to simultaneously use two of our airway visualization devices without needing separate displaying units. For example, you could use the SureSight video laryngoscope with aScope 4 Broncho or another airway visualization device.

### Streamlines work with patient imaging data

aView 2 Advance and aBox 2\* offers connectivity with Electronic Health Records (EHR) to simplify the important task of sharing and transferring patient data. You can easily view patient data before the procedure and transfer your newly recorded patient imaging to the PACS after the procedure. No additional personnel or equipment are needed. If PACS is not available in your department, you can transfer all patient imaging and procedure information files to your hospital's local network drive (SMB).

*\*VivaSight and our rhinolaryngoscopes are compatible with aView 2 Advance but not with aBox 2.*



## KEY BENEFITS

Provides clear and crisp, high contrast images

Designed for minimal set-up time, with a live image shown within 3 seconds of blade attachment

Supports dual view intubation procedures with another of our airway visualisation devices via aView 2 Advance and aBox 2

Compatible with the aView 2 Advance and aBox 2 endoscopy systems

Reliable, intuitive and high quality

Ergonomic design for secure grip and optimal control

Comprehensive range of adult blades

Designed for fast and easy blade attachment

Record images and videos at the push of a button

# Ambu

## TECHNICAL SPECIFICATIONS

### Electrical System

Power supply	Power is supplied via connection to a compatible Ambu displaying unit
User interface	1 button for image capture and video recording

### Transportation Conditions

Temperature [°C, (°F)]	-10 - 55 (14 - 131)
Relative humidity [%]	10 - 95
Atmospheric pressure [kPa, (psi)]	50 - 106 (7 - 15)

### Storage conditions

Temperature [°C, (°F)]	10 - 25 (50 - 77)
Relative humidity [%]	10 - 85
Atmospheric pressure [kPa, (psi)]	50 - 106 (7 - 15)

### Operating environment

Temperature [°C, (°F)]	10 - 40 (50 - 104)
Relative humidity [%]	30 - 85
Atmospheric pressure [kPa, (psi)]	80 - 106 (11 - 15)
Altitude [m]	≤ 2000

### General

Use (lifetime)	Reusable (3 years, or 3,120 uses)
IP Classification (SureSight system)	IP34
Biocompatible	Yes
Free of polyvinyl chloride (PVC)	Yes
Free of natural rubber latex	Yes

### Compatibility


Ambu devices	<ul style="list-style-type: none"> <li>- Ambu SureSight Blades</li> <li>- Ambu aBox 2 (Software version v.2.4.0 or newer)</li> <li>- Ambu aView 2 Advance (Software version v.2.4.0 or newer)</li> </ul>
Agents and lubricants	<ul style="list-style-type: none"> <li>- Water-based lubricants</li> <li>- Isotonic saline solution</li> <li>- 70 % ethanol</li> <li>- 70 % isopropyl alcohol</li> </ul>
Gel and solutions	<ul style="list-style-type: none"> <li>- Disodium tetraborate</li> <li>- Subtilisin</li> <li>- Isopropyl alcohol</li> <li>- Quarternary ammonium compounds</li> <li>- n-Alkyl Dimethyl Benzyl Ammonium Chloride</li> </ul>

## ORDERING SPECIFICATIONS


Item No.	Product Description	Box Quantity
390510000	Ambu SureSight Connect	1

# Ambu

 **Ambu A/S**  
 Baltorpbakken 13  
 2750 Ballerup  
 Denmark  
 T +45 72 25 20 00  
 ambu.com

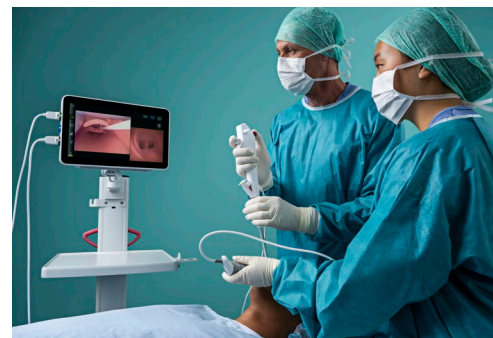

SureSight Connect can be used simultaneously with another of our airway visualisation devices for a dual view technique

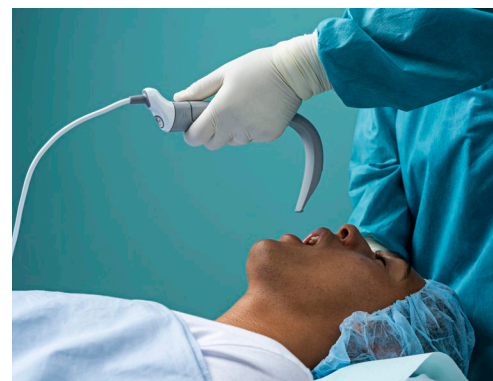

SureSight Connect has a range of compatible blades including; Macintosh, Hyper-Angulated and Channelled Hyper-Angulated options



Recording images and videos can be done at a push of a button