

DATASHEET

Ambu® Baby is a lifelike manikin that simulates babies up to the age of one year.

With the new iQF App, you can provide training with a high level of precision for infants in cardiac arrest.

CPR Training with infants is very important because it saves the lives of the youngest

Brain damage can be caused if no action is taken within four minutes without oxygen. Fortunately, this does not happen frequently. Therefore being educated in CPR training is vital.

iQF, Training that Works for Life

Helping health instructors to train. The iQF application includes data-driven information and accurate measurements of relevant CPR parameters. Help your students perform perfect CPR while enjoying games and competing with each other.



Ambu iQF Instructor

Realistic Ventilation and Ambu Hygienic System

Students will feel how a realistic airway opens when the head is properly moved into the sniff position. The instructor can simulate obstruction of the airway via an on/off slide. The Ambu Baby has additional features such as a realistic brachial pulse to provide precise training. Stomach ventilation can further be checked visually. Furthermore, the unique hygienic system prevents cross-contamination. Additionally, the correct compression depth and required ventilation volume can be monitored by the feedback indicator on the manikin.

Ambu® Baby iQF



KEY BENEFITS

- The ventilation can be seen and the exhalation through the nose can be seen and heard
- Head mobility and extension ensure exact positioning and airway clearance
- Ambu hygienic system eliminates the risk of cross infections
- Airway obstruction can be turned on via a switch
- The realism of the brachial pulse allows precise and complete teaching

TRAIN IQF, TRAIN FOR LIFE

Visualization in detail of the compression depth and the correct ventilation volume in your fingertips with Instructor and Trainee iQF App. Connected by Bluetooth®.

Ambu

SPECIFICATIONS

Weight	1,6 kg
Length	59 cm
Width	25 cm
Height	13 cm

MATERIALS

Head	Polyvinylchlorid, hard
Body skin	Polyvinylchlorid, soft
Face piece	Polyvinylchlorid, soft
Arms and legs	Polyvinylchlorid, soft
Head bag, disposable	Polyethylene
Cloth	Cotton

ORDERING SPECIFICATIONS


Item no.

A256003000	Ambu Baby iQF with 5 face pieces, 100 head bags, panties and transport bag
256000703	Ambu face piece (set of 5)
256000702	Ambu head bags, disposable (set of 100)


* Bluetooth is a registered trademark of Bluetooth SIG, inc.



Unique Ambu Hygienic System

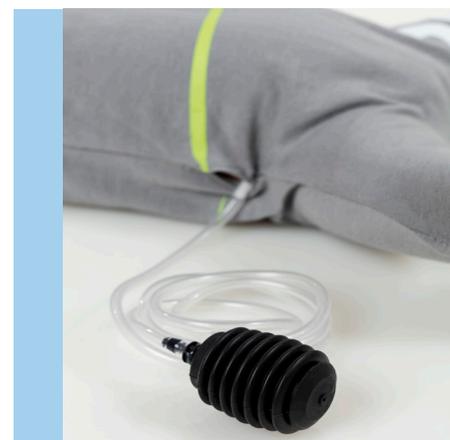

Pulse Bellow

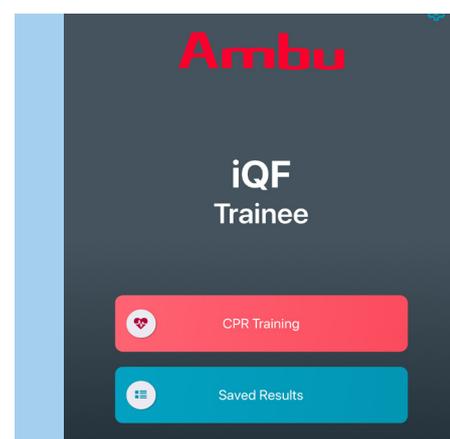

Ambu iQF Trainee

Ambu

■ Ambu A/S
Baltorpbakken 13
2750 Ballerup
Denmark
T +45 72 25 20 00
F +45 72 25 20 50
ambu.com